

DalimsurgNET

a member of the
Frankenman Group 

THE NEXT GENERATION SMART TROCAR TRANSPORT

Single-Use Trocar
Einmal Trokar Trocart
à usage unique
Trocar Monouso
Trocar de uso único
Trocares desechables



Brand Identity

TRANSPORT was developed to improve on and compensate for shortcomings that were found in existing trocars within the market. Sizes having been adjusted, now vary between 3mm to 12mm. It also has a grasp anchoring system and smoke filter(Excluding some models).



Star Sealing System

The Star Sealing System is a unique universal sealing system shaped like a star with several layered leaves to hold instruments firmly and softly.



D-Filter

TRANSPORT comes with two smoke filtering systems that filter out hazardous surgical smoke, while providing a clear vision during surgery.



T-Grip

The T-Grip was ergonomically designed as a grasp handling system for improved grip and slip prevention, in order to efficiently insert the instrument into the abdomen.



Grasp Anchoring System

The ergonomically designed Grasp Anchoring System provides better trocar fixation onto the abdominal wall, preventing accidental misplacement during surgery.



Glide Design(Glide Housing)

The Glide Housing allows the surgeon to softly glide other instruments into the curved funnel without causing damage to tissue.



Condensing Light Stinger^(Clearly Optical Viewer)



The Condensing Light Stinger provides excellent visualization and illuminates microscopic objects to provide superior contrast and detailed recognition.




D-FILTER

The smoke filtering system filters out hazardous surgical smoke, while providing a clear vision during surgery.

Close

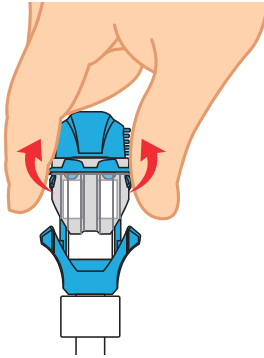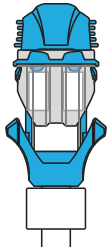


Open

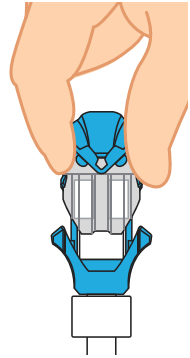

Anchoring System

The ergonomically designed Grasp Anchoring System provides better trocar fixation onto the abdominal wall, preventing accidental misplacement during surgery.

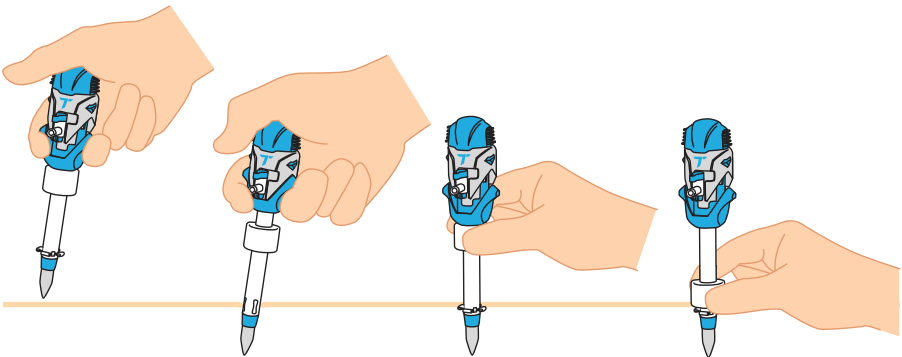

TR05F, TR05FL, TR05FL110

TR05 series is made of high quality plastic and durable materials and designed to support Reduced Port Surgery, Minimally Invasive Surgery. 5mm TRANSPORT has 3 different types of models which could be selected responding properly to the status of patients and surgeries.



Feature & Benefit

- ▶ Multi-layered silicon leaves shaped like a star to prevent air leakage.
- ▶ One smoke filtering system that filters out hazardous surgical smoke.
- ▶ Grasp handling system for improved grip and slip prevention.
- ▶ Trocar's erect fixation onto the abdominal wall. (TR05F, TR05FL110)
- ▶ Soft glide of instruments into the curved funnel.
- ▶ Visualization of tissue layers during insertion.

Model Type	Length	Diameter	Note	Quantity
TR05F	80mm	5mm	Filter, Anchor	5 units /box
TR05FL	80mm	5mm	Filter	5 units /box
TR05FL110	110mm	5mm	Filter, Anchor	5 units /box

REY MEDICAL SA: YOUR PARTNER IN SURGERY

• **QUALITY** • **SECURITY** • **ECONOMY**

Visit our website: www.reymedical.ch or call us: 022 743 24 10